

Effitec

Order them online!

Effitec OSTEOCENTRAL INFILTRATION
 100 sterile, stainless steel, siliconized, single use, dental needles with bevel mark for osteocentral and infiltration injections
 ø 0,30 (30G) 16 mm (short)

Effitec DENTALHITEC
 PATENTED SPECIAL BEVEL
 100 sterile, stainless steel, siliconized, single use, dental needles with bevel mark for osteocentral and intraligamentary injections
 ø 0,30 (30G) 9 mm (extra short)

Effitec TRANSCORTICAL OSTEOCENTRAL
 PATENTED SPECIAL BEVEL
 100 sterile, stainless steel, siliconized, single use, dental needles with bevel mark for transcortical and osteocentral injections
 ø 0,40 (27G) 16 mm (short)

Effitec NERVE BLOCK INJECTION
 PATENTED SPECIAL BEVEL
 100 sterile, stainless steel, siliconized, single use, dental needles with bevel mark for Nerve block injections
 ø 0,40 (27G) 35 mm


* Box of 100 sterile single-use dental needles

"Even previously uncomfortable injections like palatal infiltrations or upper medial injections are completely pain-free for our patients - why would I use anything else?"

Effitec needles are unique and specially designed for completely pain-free penetration into soft

tissues. There is no contact with the periosteum thanks to a very slight angle. The penetration is effective and does not cause trauma to the bone. There are four needles in the range to cover all needs."

Dr. A.PRYNNE - Bournemouth (UK)

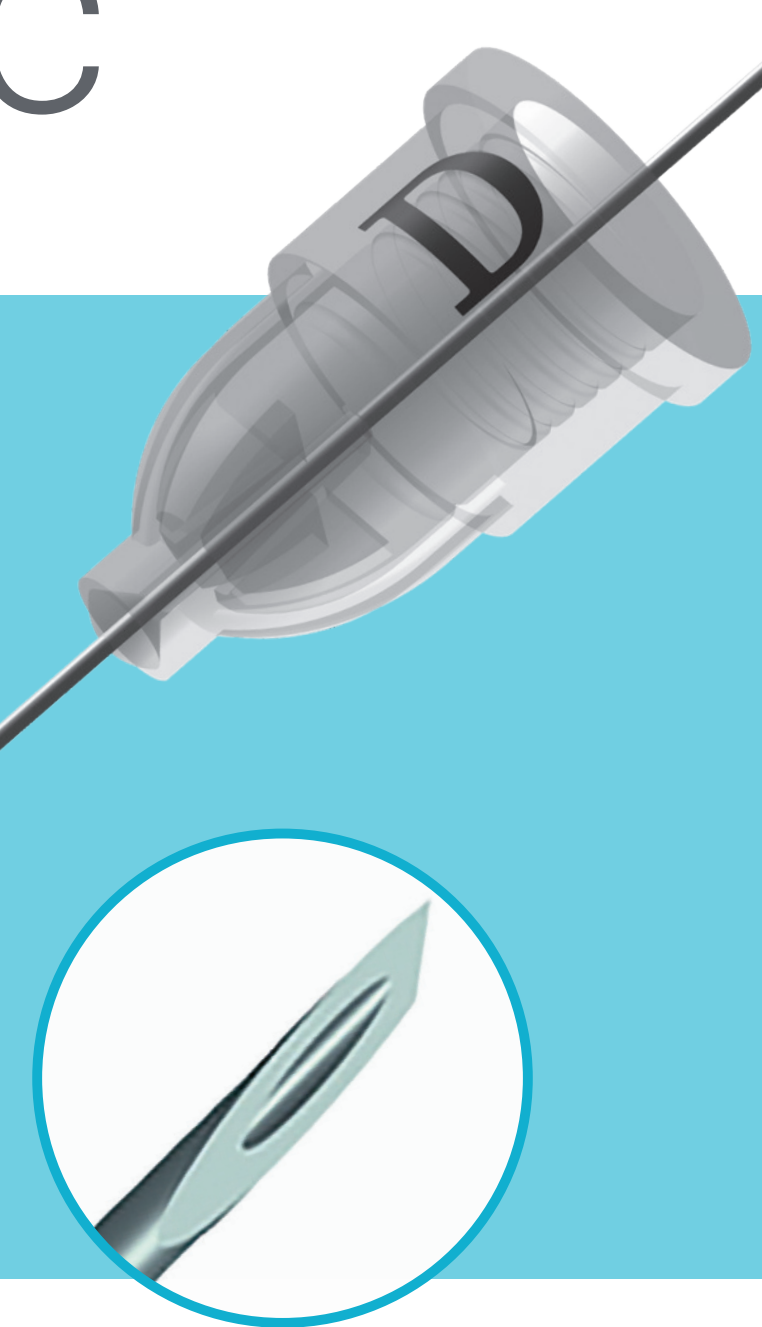

A unique scalpel bevel for guaranteed effectiveness and comfort.



ZI de l'Appentière · 49280 Mazières-en-Mauges · FRANCE
 +33 (0)2 41 56 41 91 · mail@dentalhitec.com



www.dentalhitec.com






Effitec

Why choose Effitec needles?

Their double scalpel bevel incises tissues without tearing them.

- + The needle **penetrates deeper** and without twisting for intraligamentary and intraseptal injections.
- + For attached mucosa and palatal injections, a specific and "intelligent" indicator on the hub helps to correctly orient the bevel for pain-free and minimal penetration depth.
- + Your anaesthesia will be **easier, more effective, and quicker**.
- + Needles specifically developed **for osteocentral anaesthesia** thanks to their flat bevel surface.

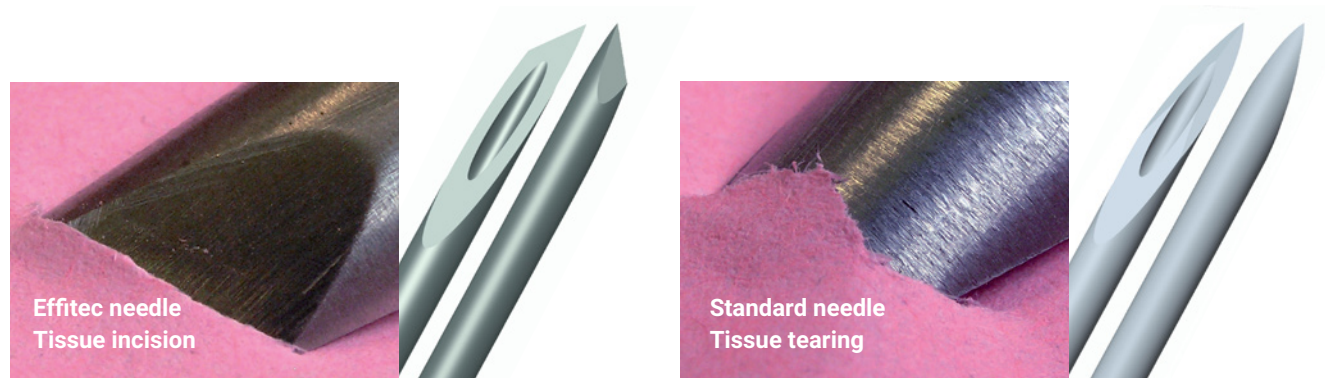


The needle revolution

"Depending on what it is used for, the needle plays a role in whether or not the injection is painful and in the success of the anaesthesia. For all of the injections made, the practitioner is first confronted with the pain caused by inserting the needle. Classic needles (single bevel) penetrate the tissues by pulling them apart and thereby create tension that generates pain. The latest generations of needles (double bevel) cut the tissue like a scalpel.

They release tension and are therefore less painful."

Dr. Alain VILLETTE

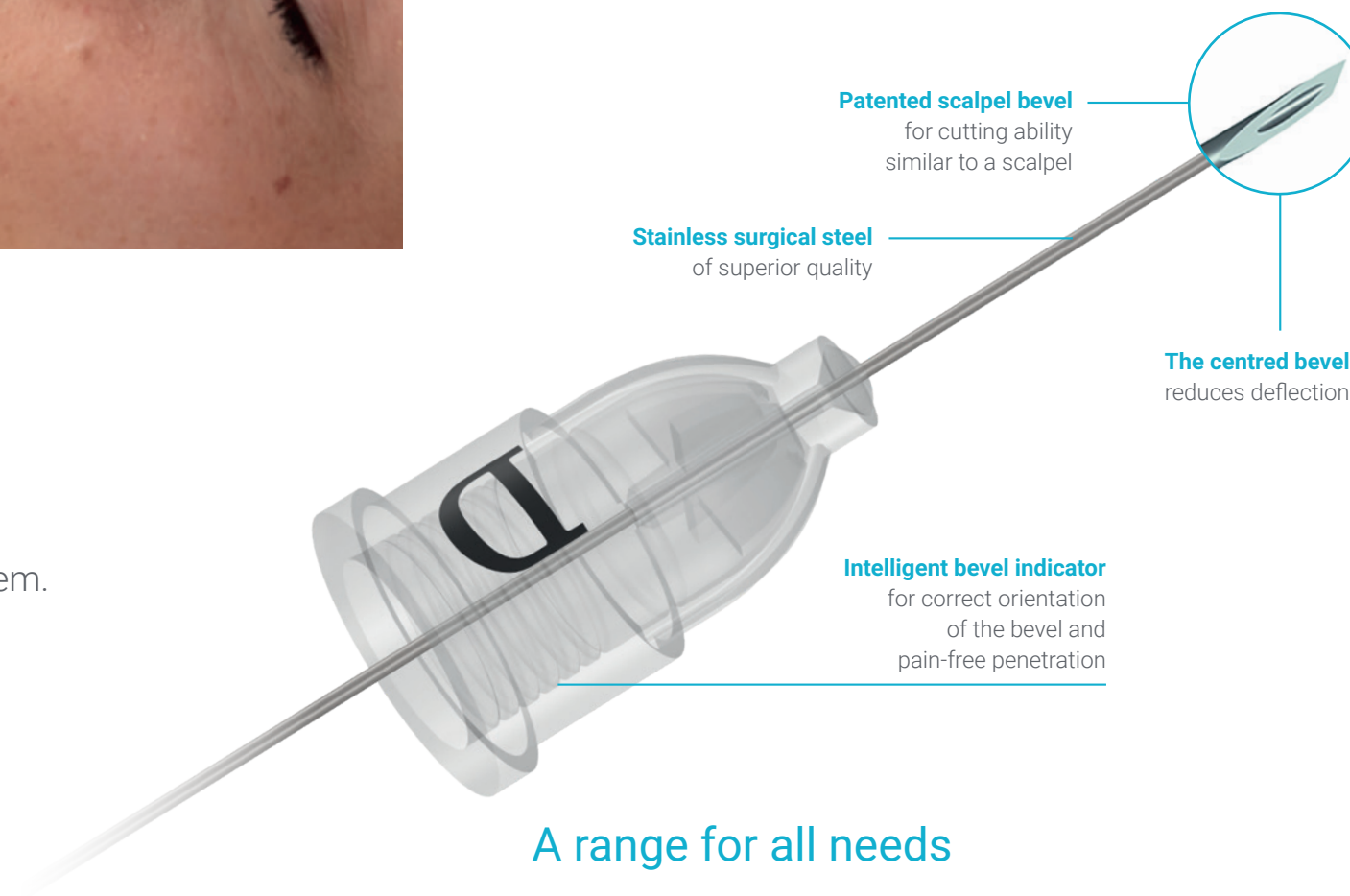

A range for all needs

Type of needle: Hub:	30G - 16 mm White	30G - 9 mm Green	27G - 16 mm Yellow	27G - 35 mm Grey
Supraperiosteal	★			
Infiltration (attached gingiva)	★	★		
Intraseptal (adults)		★		
Osteocentral	★	★ (children)	★	
Intraligamentary		★		
Transcortical			★	
Palatal	★	★		
Inferior alveolar nerve block				★

More information on www.dentalhitec.com



Easy and pain-free penetration

The mucosa is incised moving less tissue for greater accuracy and comfort.



Adapted to intrasosseous anaesthesia

The perforation is easier to perform and does not heat the bone



Compatibility

Effitec needles are compatible with all syringes on the market and with all anaesthesia techniques.