

EGO-FILLER2

WATER PURIFY SYSTEM



Provide qualified purified water source for dental equipment


- Efficient, clean and safe
- Exquisite appearance, simple design




Applicable to a variety of installation situation, beautiful and generous



Cleaning room installation diagram

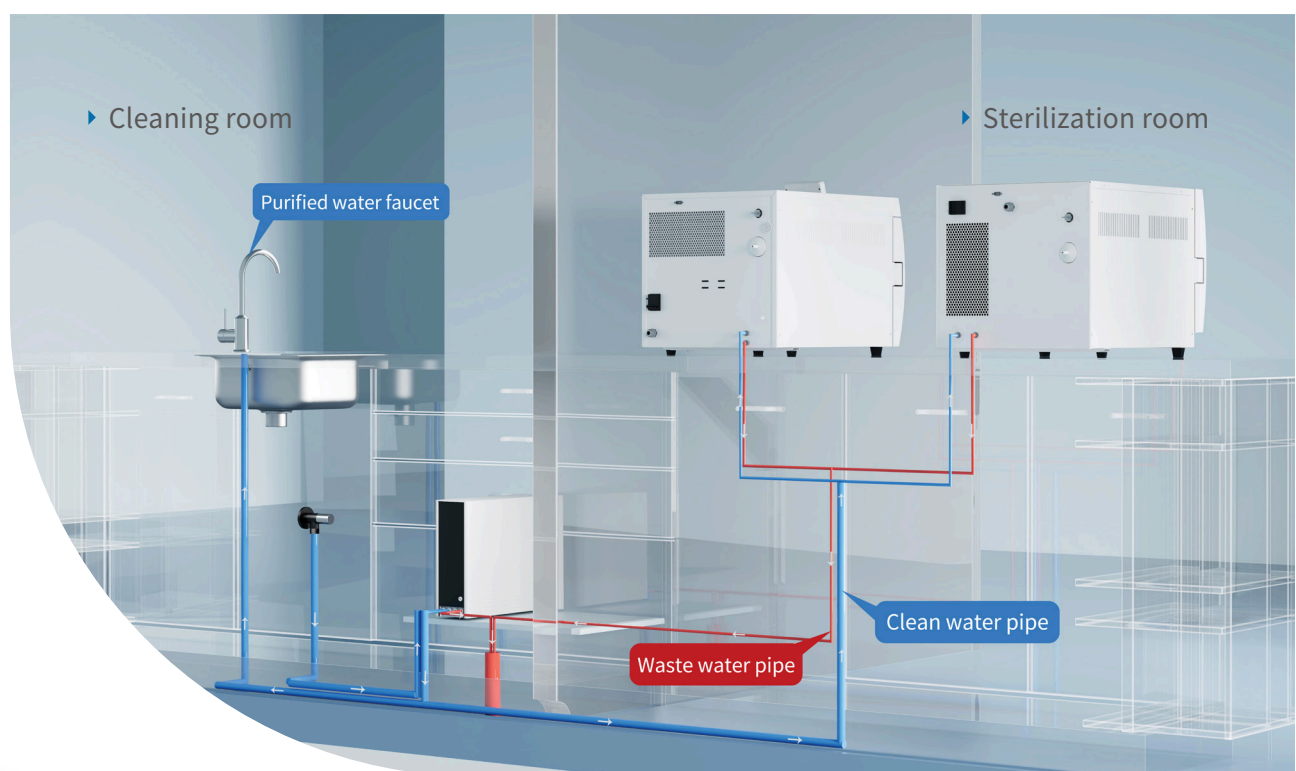


Intelligent and efficient

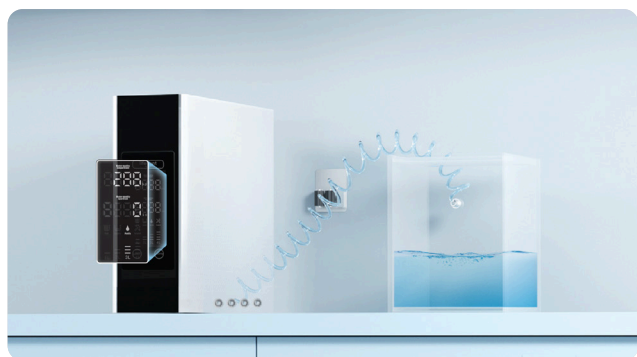

Dual water supply solution


- Freely switching automatic/manual water supply modes
- Connecting with Lifedent Apollo sterilizer to achieve automatic water supply
- Support for quantitative water supply to the water tank separately
- 1L~4L four-level water volume is selectable

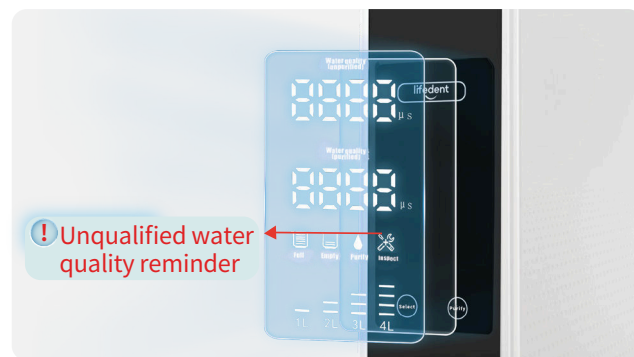

High-end digital display

- Highly sensitive touch screen, one-key water filtration
- Clear and large screen with water quality information
- Unqualified water quality reminder

Long service life Low maintenance cost

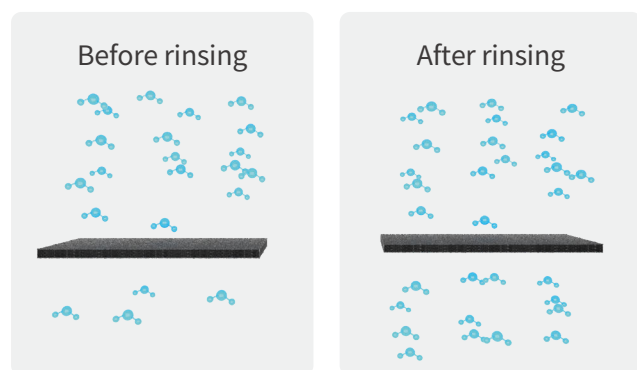

Before rinsing: Impurities will clog the surface of the filter hole and reduce the water yield

After rinsing: The surface of the filter hole is not blocked, and the water output is kept normally.

Automatic rinsing system

- Extend the service life of RO membrane and reduce the frequency of filter replacement




- Flow rate:15L/H
- It can meet the water demand of multiple sterilizers at the same time
- Total amount of purified water: 4000L(The water quality varies from place to place. Please refer to the actual situation.)
- Long service life
- Easy installation and low maintenance cost

Deep purification Clean and reliable

RO reverse osmosis


- Using high-precision RO reverse osmosis filter element, stable performance and reliable quality.
- High-efficiency RO membrane, no damage to the filter element and more water-saving.
- The filter element adopts quick connector design which is easy to be replaced.
- Built-in filter replacement reminder.

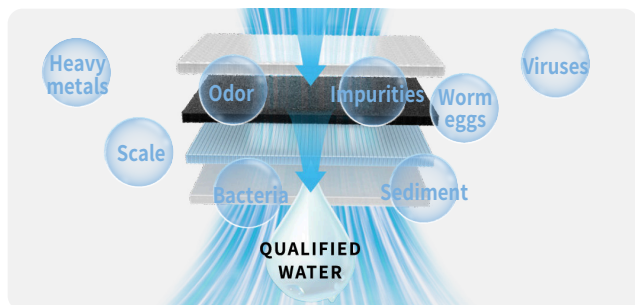


High-precision RO reverse osmosis filter element

Four-stage purification

- 4 Level filtering device which effectively filters out the rust, heavy metals, organic matter, etc.
- The deep purified water can meet the water quality standards of disinfection/sterilization.
- Avoid damage to the equipment if the water quality is not up to standard.

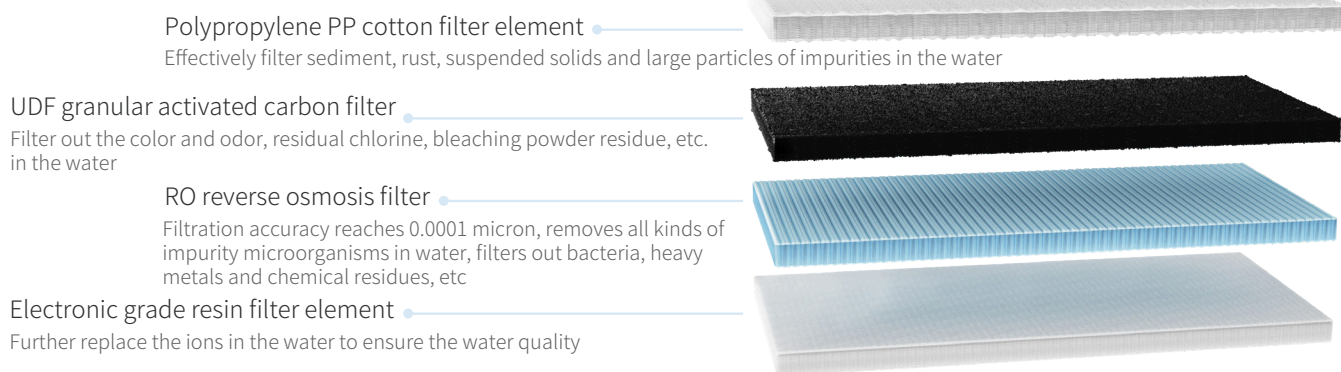

Configuration list

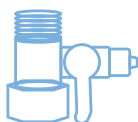

Main machine*1 Power Adapter*1 Water inlet quick connection tee*1 Spring tube*1 Counterweight connector*1 Connector holder*1 PE pipe*5m

Technical specification

Size	450*136*370mm (L *W* H)	Flow rate	15L/h
Filtration mode	Ro reverse osmosis filtration	Outlet water conductivity	<15 μ s
Filter level	Four filtering levels	Rated total amount of purified water	4000L
Display	LCD display	Applicable water quality	Municipal tap water

*The water quality varies from place to place. Please refer to the actual situation.



Lifedent(Zhuhai)Co., Ltd

Add: Floor 2 - 4, Building 1 of District A, 3588 Zhuhai Avenue, Xiangzhou District, Zhuhai, China

Tel: 86-756-8699590 Email: sales@lifedent.cn Web: www.lifedent.cn

